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Product description:

- 1 Venster en ventilatiegaten
Fenêtre
Window
- 2 Deksel
Couvercle de l'appareil
Equipment lid
- 3 Batterijencompartiment aan de onderzijde van het apparaat
Casier à piles dans le socle du boîtier
Battery compartment at the floor of the equipment
- 4 LED-indicator
Voyant D.E.L
LED indicator
- 5 UVC lichtbron
Source de lumière U.V.-C
UVC light source
- 6 Kamer voor de UVC ontsmetting en het drogen
Chambre interne pour la stérilisation U.V.-C. et le séchage
Interior chamber for UVC disinfection and drying
- 7 Startknop voor de 5 minuten ontsmetting
Bouton « ON » pour 5 min. de stérilisation
ON Start button for 5 minute disinfection

Inhoud verpakking:

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GB Cedis SanDry Eco

The equipment for the gentle disinfection and drying of the entire hearing system, preserving the integrity of the hearing system.

Drying, protecting and disinfecting in one. Your hearing system is in good hands in a Cedis SanDry Eco. It keeps your hearing system functional over time and ensures a long operating life.

The regular exposure of your hearing system with UVC light offers the utmost protection against viruses and bacteria – without liquids and additives. Even multi-resistant germs are deactivated through the UVC technology and thus become harmless. Moisture in the hearing aid, for example due to perspiration, will lead to premature deterioration of the high-quality electronics. The Cedis SanDry Eco initially disinfects and subsequently gently and safely dries the hearing aid and the earmold. This significantly increases the life span of your hearing aid.

Performance characteristics:

UVC light tube: 1.0 W
UVC wavelength: 253.7 nm (ozone free)
Power supply: DC 4.5 V – 0.35 A (3* LR 03)
LED display: blue LED-diode
Outer dimensions: W 150 x D 72 x H 40.5 mm
Inner dimensions of interior chamber:
W 85 x D 53 x H 19 mm

Operating instructions:

Use the Cedis SanDry Eco after the daily cleaning of your hearing system with Cedis cleansing wipes or cleansing spray.

1. Place the supplied 3 batteries into the battery compartment at the underside and ensure the correct poling of the batteries.
2. Open the lid of the equipment.
3. Remove the batteries from your hearing aid prior to the drying. Place the hearing system in the interior chamber with a drying capsule.
Hint: Leave the battery chamber of the hearing aid open for effective drying.
4. Press the "ON" button.
5. Close the lid of the equipment. The UVC disinfection commences for 5 minutes, the disinfection modus is signaled by the blue LED. Rapid flashing of the LED (2 times per second) indicates that the batteries should be replaced immediately to guarantee the operational readiness of the device. Upon conclusion of the UVC disinfection the LED switches off and the disinfection is completed.
6. The manual drying should occur for a further 8 hours. After opening the lid you may remove your hearing system. The capsule can be re-used until the golden-yellow color indicator has turned whitish.

Attention: Do not open capsule! Do not swallow capsule!

Hint: After every use please clean the interior chamber with a Cedis cleansing wipe.

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Important hints:

- Do not use solvents or other aggressive cleansing aids for the cleaning of your equipment.
- Do not clean the equipment under running water.
- At any sign of damage or after inappropriate use please take the Cedis SanDry Eco to your dealer for a check-up.
- The appliance may only be opened by qualified persons.
- The appliance is not to be used by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction. Children must be supervised to ensure that they do not play with the appliance.
- Do not look into the light source. Protect the skin from exposure to radiation.
- Use only for the drying of hearing aids and earmolds.
- Please dispose of batteries appropriately.
- This equipment is exclusively designed for personal use and not for medical purposes.

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Attention: In case of inappropriate usage of the Cedis SanDry Eco all warranty and liability claims are excluded.

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Operating instructions for
Cedis SanDry Premium on reverse side.


egger Otoplastik + Labortechnik GmbH
Aybühlweg 59
87439 Kempten/Germany
Tel.: +49 (0)831 58113-60
Fax: +49 (0)831 58113-14
E-Mail: sales@egger-labor.de
Internet: www.egger-labor.com

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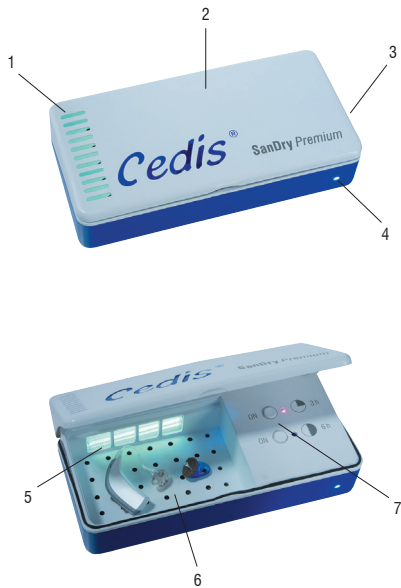Scope of delivery:
Cedis SanDry Premium
Power supply with USB connection cable and storage bag
Operating instructions

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GB Cedis SanDry Premium

The professional equipment for the gentle disinfection and drying of the entire hearing system, preserving the integrity of the hearing system.

Drying, protecting and disinfecting in one. Your hearing system is in good hands in a *Cedis SanDry Premium*. It keeps your hearing system functional over time and ensures a long operating life.

The regular exposure of your hearing system to UVC light offers optimal protection against viruses and bacteria – without liquid and additives. Even multi-resistant germs are deactivated by the UVC technology and thus become harmless. Moisture in the hearing aid, for example due to perspiration, will lead to premature deterioration of the high-quality electronics. The *Cedis SanDry Premium* initially disinfects and subsequently gently and safely dries the hearing aid and the earmold. This significantly increases the life span of your hearing aid.

Performance characteristics:
UVC light tube: 1,0 W
UVC wavelength: 253.7 nm (ozone free)
Power supply: Input AC 100–240 V ~ 50/60 Hz
Output DC 5 V === 0.5 A
Temperature range: 45 °C +/-5 °C (113 °F +/-9 °F)
LED display: LED-diode with color-change function (blue-red)
Outer dimensions: W 150 x D 72 x H 40.5 mm
Inner dimensions of interior chamber: W 85 x D 53 x H 17 mm

Operating instructions:
Use the *Cedis SanDry Premium* after the daily cleaning of your hearing system with *Cedis* cleansing wipes or cleansing spray.

1. Plug the Mini USB plug of the supplied 5 V power supply into the right-hand socket on the side of the equipment. Now plug the power supply into a suitable power socket and your *Cedis SanDry Premium* is operational. With the application of the USB cable the drying case can also be operated with suitable equipment (PC/Laptop) with USB connections.
2. Open the lid of the equipment.
3. Remove the batteries from your hearing aid prior to the drying. Place the hearing system in the interior chamber. **Hint: Leave the battery chamber of the hearing aid open for effective drying.**
4. By pressing the “ON” button 3 h or 6 h please choose the drying time for your hearing system. The chosen time will be confirmed with the lighting of the red LED diode. The lighting of the blue LED diode at the exterior of the equipment shows the operational readiness. **Hint: In case of high humidity or if you were perspiring extremely we recommend the 6 hour intensive drying.**
5. Close the lid of the equipment. The UVC disinfection commences for 5 minutes. Upon conclusion of the UVC disinfection the *Cedis SanDry Premium* automatically changes to the drying modus, which is indicated by the red LED. Upon conclusion of the drying the equipment switches off automatically; the red LED switches off. After opening the lid you may now remove your hearing system.

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6. Please allow your hearing system to cool down to body temperature for approx. 3–5 minutes after drying in the *Cedis SanDry Premium*.

Hint: After every use please clean the interior chamber with a *Cedis* cleansing wipe.

The equipment has a safety switch-off during the UVC disinfection if the lid of the equipment is opened during operation. After re-closing of the lid the UVC disinfection commences anew.

Important hints:

- Do not use solvents or other aggressive cleansing aids for the cleaning of your equipment.
- Do not clean the equipment under running water.
- Do not remove the power supply from the power plug by pulling on the cable.
- In case of damaged cable the immediate exchange is necessary.
- At any sign of damage or after inappropriate use please take the *Cedis SanDry Premium* to your dealer for a check-up.
- The appliance may only be opened by qualified persons.
- The appliance is not to be used by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction. Children must be supervised to ensure that they do not play with the appliance.
- Do not look into the light source. Protect the skin from exposure to radiation.
- Use only for the drying of hearing aids and earmolds.
- The equipment emits heat downward.
- Do not operate on heat-sensitive surfaces.
- Do not cover the equipment.
- This equipment is exclusively designed for personal use and not for medical purposes.

Warranty: 1 year from date of purchase with appropriate operation. Please submit purchasing docket.

Attention: In case of inappropriate usage of the *Cedis SanDry Premium* all warranty and liability claims are excluded.

Customer service: In case of damages or default please contact your dealer.

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Operating instructions for *Cedis SanDry Eco* on reverse side.

egger Otoplastik + Labortechnik GmbH
Aybühlweg 59
87439 Kempten/Germany
Tel.: +49 (0)831 58113-60
Fax: +49 (0)831 58113-14
E-Mail: sales@egger-labor.de
Internet: www.egger-labor.com